



## COVIDSeq as Laboratory Developed Test (LDT) for Diagnosis of SARS-CoV-2 Variants of Concern (VOC)

Rob E Carpenter<sup>1,2,5</sup>, Vaibhav K Tamrakar<sup>3,4</sup>, Sadia Almas<sup>1</sup>, Emily Brown<sup>1</sup>, Rahul Sharma<sup>1,4,5\*</sup>

### Abstract

Rapid classification and detection of SARS-CoV-2 variants have been critical in comprehending the virus's transmission dynamics. Clinical manifestation of the infection is influenced by comorbidities such as age, immune status, diabetes, and the infecting variant. Thus, clinical management may differ for new variants. For example, some monoclonal antibody treatments are variant-specific. Yet, a U.S. Food and Drug Administration (FDA)-approved test for detecting the SARS-CoV-2 variant is unavailable. A laboratory-developed test (LDT) remains a viable option for reporting the infecting variant for clinical intervention or epidemiological purposes. Accordingly, we have validated the Illumina COVIDSeq assay as an LDT according to the guidelines prescribed by the College of American Pathologists (CAP) and *Clinical Laboratory Improvement Amendments* (CLIA). The limit of detection (LOD) of this test is Ct<30 (~15 viral copies) and >200X genomic coverage, and the test is 100% specific in the detection of existing variants. The test demonstrated 100% precision in inter-day, intra-day, and intra-laboratory reproducibility studies. It is also 100% accurate, defined by reference strain testing and split sample testing with other CLIA laboratories. Advanta Genetics LDT COVIDSeq has been reviewed by CAP inspectors and is under review by FDA for Emergency Use Authorization.

**Keywords:** CLIA; Laboratory Developed Test (LDT); Next Generation Sequencing; SARS-CoV-2; Variants of Concern

### Introduction

Severe Acute Respiratory Syndrome (SARS) SARS-CoV-2 is the etiological agent of COVID-19, which is associated with mild respiratory symptoms in most infections. However, for patients with underlying medical conditions, comorbidities, and advanced age, COVID-19 may lead to severe illness. The primary route of SARS-CoV-2 transmission between humans is the respiratory route, including droplets of saliva or discharge from infected patients. Diagnosis of COVID-19 relies on detecting SARS-CoV-2 viral RNA from a nasopharyngeal or oropharyngeal specimen [1]. However, the rapid emergence of several variants with higher virulence and infectivity has provoked repeat waves of the deadly pandemic in many countries and raised anxieties about vaccine efficacy and diagnostic accuracy [2]. Rapid classification and tracking of emerging variants are critical for understanding the transmission dynamics of this disease and developing strategies for severing the transmission chain. Next-Generation Sequencing (NGS) remains the tool of choice for whole-genome analysis and deciphering new mutations [3]. Within a relatively short period, SARS-CoV-2 has acquired several mutations resulting in different virus variants. In December 2020, the United Kingdom reported a SARS-CoV-2 variant of concern (VOC),

### Affiliation:

<sup>1</sup>Advanta Genetics, 10935 CR 159, Tyler, Texas 75703, USA

<sup>2</sup>University of Texas at Tyler, 3900 University Boulevard, Tyler, Texas 75799, USA

<sup>3</sup>ICMR-National Institute of Research in Tribal Health, Jabalpur, MP 482003, INDIA

<sup>4</sup>RetroBioTech LLC, 838 Dalmalley Ln, Coppell, TX 75019, USA

<sup>5</sup>Scienetix, 10935 CR 159, Tyler, Texas 75703, USA

### \*Corresponding author:

Rahul Sharma, Advanta Genetics, 10935 CR 159 Tyler, Texas 75703, USA.

**Citation:** Rob E Carpenter, Vaibhav K Tamrakar, Sadia Almas, Emily Brown, Rahul Sharma. COVIDSeq as Laboratory Developed Test (LDT) for Diagnosis of SARS-CoV-2 Variants of Concern (VOC). *Archives of Clinical and Biomedical Research* 6 (2022): 954-970.

**Received:** November 11, 2022

**Accepted:** November 21, 2022

**Published:** December 02, 2022

lineage B.1.1.7, detected in over 30 countries and is more efficiently transmitted than other SARS-CoV-2 variants. Thus, the pandemic strikes in several phases of outbreaks in different parts of the world [4]. Currently, the virus continues to be a global agent of infection. The highly mutagenic nature of SARS-CoV-2 assaulted many countries with second or third waves of the outbreak [5, 6]. Mutations with higher transmissibility, a more intense disease state, and less likely to respond to vaccines or treatments have been classified by the World Health Organization (WHO) as Variants of Concern. Recent epidemiological reports released by WHO indicated five VOCs: 1) B.1.1.7 (Alpha) in December 2020; 2) B.1.351 (Beta) in December 2020; 3) P.1 (Gamma) in January 2021; 4) B.1.617.2 (Delta) in December 2020, and 5) B.1.1.529 (Omicron) (<https://www.who.int/en/activities/tracking-SARS-CoV-2-variants/>). The receptor-binding domain (RBD) of coronavirus increases its capacity to strike in several outbreak phases in different parts of the world [7]. More recently, South Africa reported a new SARS-CoV-2 variant to the WHO. Omicron (B.1.1.529) was first detected in specimens collected in Botswana and designated as the fifth VOC [8] (<https://www.who.int/en/activities/tracking-SARS-CoV-2-variants/>). Several variant-specific treatment options have been approved by the U.S. Food and Drug Administration (FDA), including *Bebtelovimab*, a monoclonal antibody for the treatment of COVID-19 that retains activity against the omicron variant. However, a recent study shows that the effectiveness of mRNA vaccines is reduced against all three subvariants of omicron [9]. Several other studies have reported a substantial decrease in neutralizing antibody titers after vaccination against all coronavirus variants. [10, 11]. Reduced neutralizing activity against the B.1.1.7 (Alpha), B.1.35 (Beta), and P.1 (Gamma) strains have been reported among the Pfizer-BioNTech vaccinated populations [12]. Another study investigated the neutralization of antibodies elicited by Novavax NVX-CoV2373, a protein subunit vaccine, and that of mRNA-1273 by Moderna against the California variant B.1.429 and B.1.351 pseudoviruses. The small-scale study of 63 volunteers similarly revealed the reduction in neutralization abilities of antibodies elicited by both vaccines. The most drastic reduction, up to a 9–14 times decrease in neutralization compared to D614G, was observed with B.1.351 pseudovirus, where the antibodies were 2–3 times less sensitive against the B.1.429 variant pseudovirus [13]. SARS-CoV-2 is likely to continue to evolve, and the next strain may have a strain-specific etiology requiring strain information for patient care. In the present situation, most infections are attributed to a single sublineage. However, new lineages are likely to emerge and replace existing circulating lineages. Several PCR-based assays are available for the detection of the known variants. These assays are not designed to detect unknown infected variants [14]. Unlike the PCR-based test, this NGS-based assay can detect new variants as they emerge. And because lineage variance has

potential implications for virulence and infectivity, validation of NGS assays that proactively identify mutagenic variants enables these test results to be used in clinical applications when warranted. Furthermore, this Illumina COVIDSeq assay is used for epidemiological surveillance globally. Still, the validation of the assay as a Laboratory Developed Test (LDT) is required to use variant information for clinical decision-making. For example, variant-specific monoclonal antibody therapies have been emphasized by the National Institutes of Health (NIH) COVID Treatment Guidelines Panel and the FDA, recommending against the use of bamlanivimab and etesevimab (administered together) and REGEN-COV (casirivimab and imdevimab) because of significantly reduced activity. Consequently, we report the validation of the NGS-based test to identify the existing and emerging variants of SARS-CoV-2 [16]. This study has benchmarked the validation process for using the variant information in clinical management as required by CLIA. Although the Illumina COVIDSeq assay has been approved for emergency use authorization (EUA) for the diagnosis of COVID-19; the assay has not been approved for variant detection. We validated the Illumina COVIDSeq assay according to CLIA/CAP requirements for LDT, and the validation report has been submitted to the FDA for EUA, and reviewed by a team of CAP inspectors. Accordingly, the COVIDSeq assay is qualified to diagnose SARS-CoV-2 variants of infected individuals and can be deployed for monitoring the evolution of SARS-CoV-2 variants in decentralized clinical laboratory settings.

## Materials and Methods

The workflow consists of the following procedures: RNA extraction, cDNA synthesis, target amplification, library preparation, library pooling, sequencing, and analysis. Validation was performed to achieve a high degree of accuracy and precision. Additional studies were performed to test the effect of interference substances and sample stability.

### Reference Strains of SARS-CoV-2 Variants

We used three reference strains of SARS-CoV-2, Omicron, Delta, and Wuhan. Complete genome synthetic RNAs of these strains were obtained from BEI Resources.

### Clinical Specimens

De-identified sample remnants from nasopharyngeal swabs were collected from the patients who tested positive for SARS-CoV-2 PCR with RT-PCR at Advanta Genetics (<https://aalabs.com/>) in Tyler, Texas. Samples were stored at -80°C until RNA extraction. The study was exempted by IRB (Institutional Review Board) because only de-identified samples were used.

### RNA Extraction

Total RNA was extracted using the Roche *MagNA*

Pure 96 System and Viral RNA Small Volume Kits per the manufacturer's (Port Scientific Inc. QC J3G 4S5 Canada) instructions. Isolated RNA was frozen at -80°C until the library preparation.

### Library Preparation and Sequencing

The libraries were prepared using the Illumina COVIDSeq protocol (Illumina Inc, USA). Briefly, total RNA was primed with random hexamers, and first-strand cDNA was synthesized using reverse transcriptase. The SARS-CoV-2 genome was amplified using the two sets of primers (COVIDSeq Primer Pool-1 & 2) provided by Illumina, but the primer sequences have not been disclosed by the manufacturers. Primers are not mutation specific but designed to amplify the entire genome. PCR amplicons were tagged using the EBLTS (Enrichment BLT), which is a process that fragments and tags the PCR amplicons with adapter sequences. Adaptor ligated amplicons were further amplified using the distinct pre-paired 10 base pair Index 1 (i7) adapters and Index 2 (i5) adapters (IDT for Illumina-PCR Indexes Set 1) for each sample. The individual library was quantified using a Qubit 2.0 fluorometer (Invitrogen, Inc.) and pooled in equimolar concentration instead of equal volume as recommended by Illumina. This additional step allowed us to achieve uniform coverage of all the libraries in the pool and efficient use of a low throughput sequencing instrument (MiniSeq®). A COVIDSeq positive control (Wuhan-Hu-1) and one no template control (NTC) were processed with each batch of libraries. The final library pool was again quantified using a Qubit 2.0 fluorometer (Invitrogen, Inc.) and a PCR-based library quantification kit (Scienetix, USA). The final library pool was diluted to a 2 pM loading concentration. Dual indexed paired-end sequencing with 75bp read length was carried out using the HO flow cell (150 cycles) on the Illumina MiniSeq® instrument.

### NGS Data Analysis

Illumina Basespace (<https://basespace.illumina.com>) bioinformatics pipeline was used for sequencing QC, FASTQ generation, genome assembly, and identification of SARS-CoV-2 variants. Briefly, the raw FASTQ files were trimmed and checked for quality (Q>30) using the FASTQ-QC application within the Basespace. QC passed FASTQ files were aligned against the SARS-CoV-2 reference genome (NCBI Reference Sequence NC\_045512.2) using Bio-IT Processor (Version: 0x04261818). Then, DRAGEN COVID Lineage (Version: 3.5.4) application in Bbasespace was used for SARS-CoV-2 variant determination and generating a single consensus FASTA file. Finally, single consensus FASTA was also analyzed for lineage assignment using the web version of Phylogenetic Assignment of Named Global Outbreak Lineages (PANGOLIN) software (<https://pangolin.cog-uk.io>). Only the consensus variants identified by both applications were used for further analysis.

### Strain Typing of SARS-CoV-2 in the East Texas Region

We have applied the SARS-CoV-2 variant detection workflow established in this study for strain typing of SARS-CoV-2 in the East Texas USA region over the course of the pandemic. Representative samples collected at various po time during the pandemic (Aug 2020, July 2021, Dec 2021, April 2022, July 2022, and Sept 2022) were sequenced and analyzed for the circulating variants in the region of interest.

### Results

The SARS-CoV-2 sequencing test was validated as LDT according to the guidelines prescribed by the CAP and mandated by CLIA. Briefly, the limit of detection (LOD), analytical accuracy, precision, and sample stability was established. The effect of carryover and interference substances was also investigated.

#### Analytical Sensitivity/Limit of Detection

This assay is not intended to diagnose the SARS-CoV-2 infection but is meant for discovering the SARS-CoV-2 variant from a patient previously diagnosed with a high level (Ct<30) of SARS-CoV-2 infection. LOD of this assay was determined for two variables needed for accurate results: 1) the lowest amount of input genomic material and 2) minimum genomic coverage. For genomic material LOD, serial dilutions of an Omicron reference variant were sequenced in triplicates, and the lowest input concentration resulting in the correct variant detection is identified as LOD (Table-1). LOD for the test is defined as Ct<30 (~15 copies/ul of RNA). LOD was further verified by sequencing 23 samples with RNA input close to LOD (PCR Ct value ~30±2) and obtained 200X-1000X coverage; variants for all the 23 samples were identified correctly. We have also analyzed the 26 additional samples from patients found positive during September 2022—all the samples have Ct<30. We sequenced the 26 samples, including 10-fold and 100-fold dilutions of three representative samples. We were able to identify the variant in all 26 samples and endorse the application of the assay on current circulating strains post-vaccination. However, the application of this test in asymptomatic or very low viral load (<150 viral genome/ul) remains the limitation of the assay.

To determine the LOD regarding genomic coverage, we computed the depth of coverage (X times) and percent genome coverage for all tested samples. The lowest genomic coverage of >200X (Depth) and 90% genome coverage is required for successful detection variant detection (Figure-1). Importantly, all 164/164 (100%) observations with a minimum of 90% genome coverage at a minimum of 200X resulted in the correct variant call after the analysis.

#### Analytical Accuracy

Accuracy is a determination of the amount of systematic error in the system. The analytical specificity of this assay is

determined by re-sequencing the already sequenced reference strains of the SARS-CoV-2 virus and the alignment of the resulting FASTQ files to the available reference genome sequence using the BaseSpace (Illumina) tool. We tested 3 known SARS-CoV-2 variants; Wuhan-Hu-1, B.1.617.2 (Delta), and B.1.1.529 (Omicron) in triplicates, and all the variants were identified correctly as expected. NGS does not use analyte (i.e., SARS-CoV-2 variant) specific reagents to determine the correct variant but uses whole genome analysis for discerning the variants. Therefore, the specificity of the variant detection is considered 100% (Table 2).

Considering the limited availability of reference strains, we also re-sequenced the 6 samples already sequenced by another reference laboratory (Fulgent Genetics) at extremely high coverage (>30,000X), and variant identities were compared between two observations. A total of 6 samples were sequenced at Fulgent Genetic and Advanta Genetics. All 6 samples were identified to carry identical variants by both laboratories implicating 100% accuracy in inter-laboratory testing (Table-3). The average sequencing coverage at Fulgent Genetics is 34560.5X compared to 174.3X at Advanta. Interestingly, variants of 3 samples sequenced at >50,000 coverage were correctly identified by only 200X coverage. With pre-pooling quantification, we achieved higher sequencing efficiency without compromising the test accuracy. Such higher efficiency is critical for the cost-effective application of this test in limited-resourced and de-centralized laboratory settings.


## Precision

The precision of a measurement system, related to reproducibility and repeatability, is the degree to which repeated measurements under unchanged conditions show the same results. Inter-day precision is determined by sequencing 9 samples of known genomic variants over three days. Nine samples were tested in three rounds of library preparations, sequencing, and data analysis. The identity of the variant detected across the three runs was compared. All 9 samples were identified correctly across the three sequencing instances implicating 100% precision across the 27 observations. Inter and intra-day precision was determined by testing 6 clinical samples (near LOD) in triplicate during three rounds of library preparations, sequencing, and data analysis. The identity of the variant detected across the three runs was compared. Three samples that failed the pre-defined QC (library yield, reads, coverage, etc.) were excluded from the precision. Over three days, the remaining 51/54 observations were in 100% concordance for triplicate testing. All 6 samples were identified as the same variant in triplicate testing within a single batch implicating 100% intra-run precision (Supplementary-2). Likewise, the same samples resulted in identical variants when tested in three distinct batches of library preparation, sequencing, and data analysis (Supplementary-3). Thus, inter-day precision was also determined as 100%.

**Table 1:** Limit of Detection: Four serial dilutions of Omicron strain were sequenced in duplicate, and the lowest viral RNA input, which resulted in accurate variant detection, was accepted as LOD for the test.

Pango Lineage	Sample ID	Relative Copy/ul	ACSQ Observed Lineage					
			PCR Ct Value		Median Coverage	Coverage >= 30x	Pango Lineage	WHO label
			N1	N2				
BA.1	ACSQ4-19	1.48E+04	26.16	26.37	663	96.79%	BA.1	Omicron
	ACSQ4-19 (1:10)	1.48E+03	29.87	30.15	690	95.18%	BA.1	Omicron
	ACSQ4-19 (1:100)	1.48E+02	33.24	33.54	504	94.46%	BA.1	Omicron
	ACSQ4-19 (1:1000)	1.46E+01	36.54	35.47	586	91.29%	BA.1	Omicron
BA.1	ACSQ4-19	1.48E+04	26.16	26.37	533	95.87%	BA.1	Omicron
	ACSQ4-19 (1:10)	1.48E+03	29.87	30.15	489	94.89%	BA.1	Omicron
	ACSQ4-19 (1:100)	1.48E+02	33.24	33.54	745	94.97%	BA.1	Omicron
	ACSQ4-19 (1:1000)	1.46E+01	36.54	35.47	736	93.61%	BA.1	Omicron
BA.4.6	ACSQ9-1A	1.48E+03	29.83	28.36	1597	94.25%	BA.4.6	Omicron
	ACSQ9-1A (1:10)	1.48E+02	33.16	31.69	330	82.33%	BA.4.6	Omicron
	ACSQ9-1A (1:100)	1.46E+01	36.49	35.02				Not detected
BA.5.2.1	ACSQ9-4A	1.58E+03	27.4	26.17	2300	96.89%	BA.5.2.1	Omicron
	ACSQ9-4A (1:10)	1.58E+02	30.73	29.5	1635	92.93%	BA.5.2.1	Omicron
	ACSQ9-4A (1:100)	1.58E+01	34.06	32.83	514	88.59%	BA.5.2.1	Omicron
BA.2.3	ACSQ9-23	1.48E+04	26.65	23.71	3355	99.46%	BA.2.3	Omicron
	ACSQ9-23 (1:10)	1.48E+03	29.98	27.04	2798	98.52%	BA.2.3	Omicron
	ACSQ9-23 (1:100)	1.48E+02	33.31	30.37	2051	95.68%	BA.2.3	Omicron



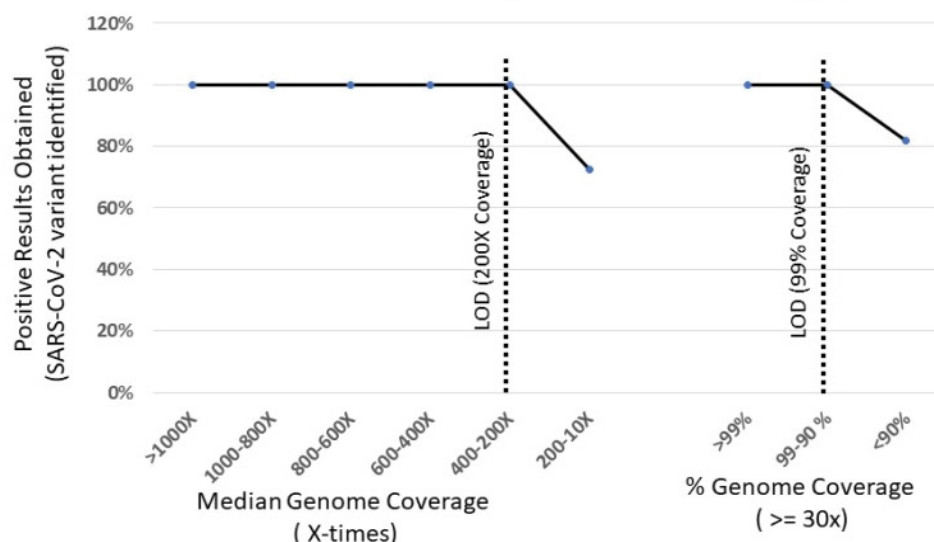

**Figure 1:** Limit of Detection (LOD): Median genomic coverage (X-times) and minimum % length of the genome covered >30X times were computed, and the minimum coverage required for obtaining the accurate SARS CoV-02 lineage was defined as LOD.

**Table 2:** Accuracy of the test is determined by sequencing 3 reference strains of the SARS-CoV-2 virus.

Sample ID	Reference Lineage			Detected Lineage			Concordance
	Median Coverage	Pango Lineage	WHO label	Median Coverage	Pango Lineage	WHO label	
Wuhan (1500cp/ul)	NA	A1	Non-VOC	1167	A1	Non-VOC	YES
Wuhan (150cp/ul)	NA	A1	Non-VOC	2289	A1	Non- VOC	YES
Omicron (1500cp/ul)	NA	B.1.1.529	Omicron	1372	B.1.1.529	Omicron	YES
Omicron (150cp/ul)	NA	B.1.1.529	Omicron	829	B.1.1.529	Omicron	YES
Delta (1500cp/ul)	NA	B.1.617.2	Delta	663	B.1.617.2	Delta	YES
Delta (150cp/ul)	NA	B.1.617.2	Delta	533	B.1.617.2	Delta	YES

**Table 3:** Comparative genome sequencing and variant calling results obtained from two different laboratories.

Sample ID	Fulgent Genetics*			Advanta Genetics			Concordance
	Median Coverage	Pango Lineage	WHO label	Median Coverage	Pango Lineage	WHO label	
ACSQ1-2	77387.8	AY.4	Delta	214	AY.25	Delta	YES
ACSQ1-3	341.6	AY.3	Delta	228	AY.3	Delta	YES
ACSQ1-4	54823.9	AY.3	Delta	222	AY.3	Delta	YES
ACSQ1-5	2860.5	AY.4	Delta	218	AY.3	Delta	YES
ACSQ1-6	69088.9	B.1.617.2	Delta	258	B.1.617.2	Delta	YES
ACSQ1-7	2860	AY.3	Delta	106	AY.3	Delta	YES

\*Variant from AY4 to AY 25 were reported as AY4 in the sequencing results obtained from the Fulgent Genetics.

## Stability Study

The stability of clinical samples at different temperatures was tested to simulate the temperature conditions during transportation. Samples identified as Omicron (n=3) and Delta (n=1) were placed at 4 different temperatures [Freezer (-20°C); Refrigerator (2-8°C); Room Temp (~25°C); Elevated

Temp (~50°C)] to mimic the possible environmental conditions during the transportation. Samples were left for up to 7 days under these temperature conditions. Samples were retrieved at 24 hours, 3 days, and 7 days intervals, and RNA was extracted and stored at -80°C. RNA from all the samples in the stability study was tested in a single library preparation

and sequencing batch. We were able to sequence and identify the SARS-CoV-2 variant of the samples kept at 20°C, 2-8°C, ~25°C, and ~50°C for 24 hours and 3 days. However, samples placed in elevated temperature conditions resulted in low-quality sequencing data, which did not result in variant detection. Overall, samples kept at an elevated temperature (~50°C) over 3 days were unsuitable for variant detection by whole-genome sequencing (WGS). All the samples used for the stability study were of viral load close to LOD.

### Freeze-Thaw Stability Study

Extracted RNA was subjected to 2 and 3 freeze-thaw cycles, and RNA was processed in single Library preparation. RNA sample after 2 freeze-thaw cycles fail the pre-defined sequencing QC and could not be used for variant detection. Although sample or RNA storage conditions are unlikely to change the SARS-CoV-2 variant, >3 days of storage at high temperature (~50°C) may cause the SARS-CoV-2 variant testing to fail or result inconclusive because of compromised data quality. Likewise, >2 freeze-thaw cycles for RNA also compromised the sequencing data quality. Thus, samples for SARS-CoV-2 variant detection should be kept at 4°C for 7 days and stored at -20°C for the long term (Table-4).

### Role of Interference Substances

Clinical specimens may contain biological or non-biological substances which may interfere with the testing process. We spiked the commonly used nasal sprays into the clinical specimen and tested the sample with and without the external substance. None of the tested substances altered the results or compromised the data quality (Data not shown).

### Epidemiological Survey of East Texas, US

We have also applied Advanta Genetics LDT COVIDSeq to investigate the evolution of SARS-CoV-2 in the East Texas region during the pandemic. We identified a greater genomic diversity in early pandemics before identifying variants of Concern. We identified the SARS-CoV-2 variant (B.1, B.1.126, B.1.2, B.1.234, B.1.243, B.1.564, B.1.574, B.1.602) among the samples collected in July 2020. All of these variants were categorized as non-VOC by the WHO. Diverse non-VOC strains were initially replaced by the Delta variant (100%) in July-Aug 2021. Omicron (58%) and Delta (42%) variants were co-circulating during Dec 2022; Delta was completely replaced by the Omicron variant by December 2021. Omicron BA.2 (79%) was the dominant variant during April 2022, which was again replaced by BA.5 (78%) in September 2022. (Figure-2). Thus, continuous monitoring is warranted to keep the pandemic from returning to the scale seen earlier by identifying the vaccine escape or target dropout in diagnostic testing. All the SARS-CoV-2 whole-genome sequences generated in this study were submitted to GISAID (<https://www.gisaid.org>) database (Supplementary Table-1), [17].

### Discussion

WHO has been classifying the SARS-CoV-2 variant into various categories according to their possible clinical implication and public health concern. Several technologies have been adopted for SARS-CoV-2 variant detection, but NGS remains the gold standard because of comprehensive genomic analysis [18, 19]. In June 2020, the US FDA granted EUA for Illumina's NGS test for COVID-19 diagnosis. However, the test has not been widely adopted for diagnosis because RT-PCR is much cheaper and easy to implement in the unprecedented need for SARS-CoV-2 testing. Although RT-PCR remains the method of choice for routine diagnosis, Illumina COVIDSeq protocol has been instrumental in outbreak investigation and surveillance throughout the pandemic [20]. Several laboratories worldwide use WGS for high throughput surveillance communicated by health organizations [21]. The Delta variant has been associated with greater transmissibility and higher viral RNA loads in both unvaccinated and fully vaccinated individuals [22]. WGS has also identified the potential compromise vaccine effectiveness against the Omicron variant [23]. The emergence of SARS-CoV-2 variants with significantly different clinical implications accentuates the need for variant detection, especially for immunocompromised patients. Additionally, because some monoclonal antibody treatments are variant-specific, timely identification of the infecting SARS-CoV-2 variant may influence decision-making and treatment. Currently, there is no FDA-approved SARS-CoV-2 variant detection test for diagnosing individual patients. Thus, LDT remains the only viable option to leverage NGS methods for SARS-CoV-2 variant diagnosis. This virus is predicted to mutate continuously, and the evolution of variants with significantly different clinical interventions cannot be ruled out [5]. This study established the WGS workflow for detecting SARS-CoV-2 variants according to CLIA guidelines for LDTs. The importance of reference materials for the validation and QC of wet-lab and dry-lab WGS processes is well established [24]. However, unlike human genomics [25], there is no well-established resource of reference materials for the validation of such genetic variant detection tests. Therefore, we obtained three reference strains (Wuhan, delta, and Omicron) of SARS CoV-2 for accuracy study, and all three variants were correctly identified in repeated testing. The Wuhan strain of SARS-CoV-2 has been accepted as the reference strain [26]. Therefore, all the sequences generated during this study were aligned against the Wuhan strain genome [26]. Although more than a million SARS-CoV-2 genomes have been sequenced [27], a limited number of well-characterized reference genomic materials are available. To overcome this limitation, we re-sequenced 6 samples already sequenced by another laboratory (Fulgent Genomics), and the sequencing results were in 100% concordance.


Interestingly, 3 of these samples were sequenced at

>50,000X coverage by Fulgent Genetics, and the same samples were sequenced at 200-300X coverage at Advanta Genetics. Results from both sequencings were in 100% concordance, suggesting that such high sequencing depth is unnecessary for routine variant detection. The introduction of pre-pooling quantification and equimolar pooling enabled us to achieve uniform distribution of sequencing reads across the samples in the pool, resulting in more efficient sequencing. This approach particularly important in de-centralized reference laboratories which do not have access to high throughput instruments such as Illumina HighSeq or NovSeq. We were able to sequence up to ~30 samples in a single MiniSeq run,

reducing the cost of sequencing (excluding library preparation) to ~ \$30/sample. The genome sequences available from public databases may have been generated using different sequencing chemistries or platforms, which might yield different error rates; therefore, the inter-laboratory study was of particular interest because the reference laboratory used a different instrument for sequencing. Overall, we achieved high accuracy, reproducibility, repeatability, diagnostic (variant detection) sensitivity, and specificity of 100%, which exceeds the 90% threshold for LDT performance parameters per CLIA requirements. These findings agree with other reports of 93% to 100% accuracy in WGS identification and

**Table 4:** Samples stored in simulated environmental conditions mimicking the possible transportation and storage temperature were sequenced to identify acceptable sample storage conditions.

Sample-ID	Temp	Time	Median Coverage	Coverage >= 30x
ACSQ4-1	Freezer (-20C)	24hr	1256.5	98%
ACSQ4-9	Freezer (-20C)	48hr	1286	98%
ACSQ5-1	Freezer (-20C)	72hr	588.5	96%
ACSQ5-9	Freezer (-20C)	7days	524	95%
ACSQ4-3	Room Temp (~25C)	24hr	1284.5	97%
ACSQ4-11	Room Temp (~25C)	48hr	1254.5	97%
ACSQ5-3	Room Temp (~25C)	72hr	643	96%
ACSQ5-11	Room Temp (~25C)	7days	458.5	95%
ACSQ4-2	Refrigerator (2-8C)	24hr	998.5	94%
ACSQ4-10	Refrigerator (2-8C)	48hr	1086.5	96%
ACSQ5-2	Refrigerator (2-8C)	72hr	537	94%
ACSQ5-10	Refrigerator (2-8C)	7days	477	94%
ACSQ4-4	Elevated Temp (~50C)	24hr	866.5	94%
ACSQ4-12	Elevated Temp (~50C)	48hr	1029	96%
ACSQ5-4	Elevated Temp (~50C)	72hr	363	91%
ACSQ5-12	Elevated Temp (~50C)	7days	Low Coverage	
ACSQ5-16	Elevated Temp (~50C)	7days	Low Coverage	

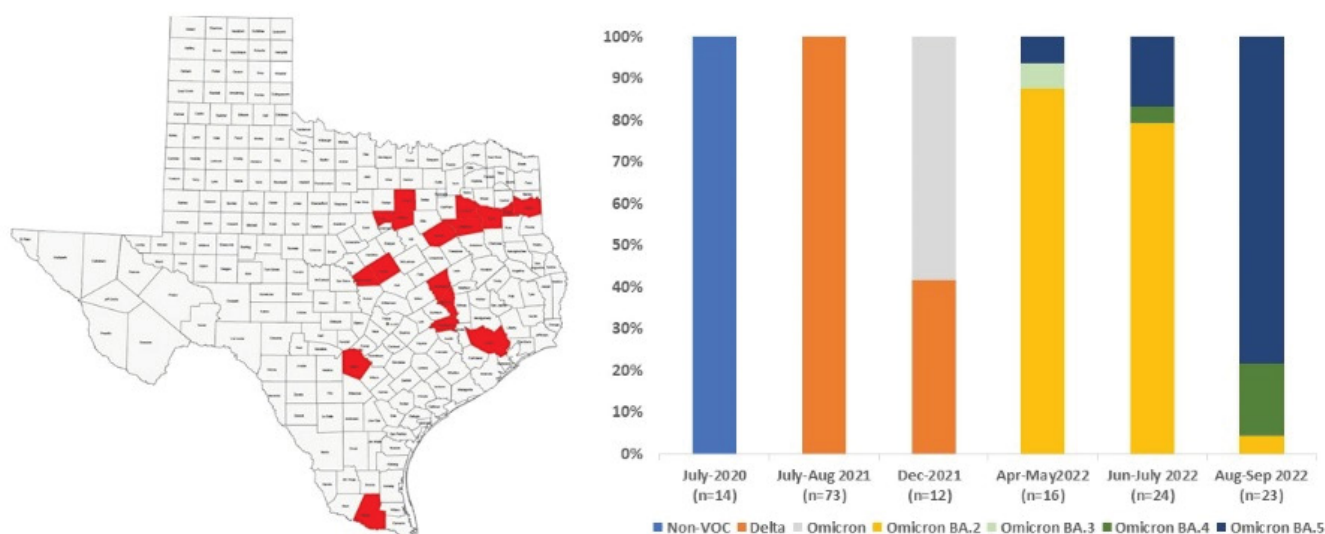

**Figure 2:** Evolution of SARS-CoV-2 variants in East Texas over the course of the Pandemic.

subtyping for other pathogens [28, 29]. We determined the LOD as 90% genome sequenced at >30X depth and >200X median depth of coverage. The LOD study did not consider coverage for individual single nucleotide polymorphisms (SNP) because an SNP combination determines the SARS-CoV-2 genomic variant. LOD, in terms of minimum genomic copies, was established at ~15 copies/μL going into the sequencing reaction. This assay can identify the genomics variant from the lower viral RNA input, but this is the lowest input tested during this validation. Vaccination has reduced hospitalization and deaths in COVID cases, but the viral load (10E+05 to 10E+08 genomic copies/ml) in breakthrough cases remains high enough for detection by this assay [15]. Interestingly, only infectious viral load (VL) was lower in fully vaccinated Omicron BA.1-infected individuals compared to vaccinated Delta-infected individuals, indicating variant-specific response to the vaccines. A reduced infectious VL was observed only in boosted but not fully vaccinated individuals compared to unvaccinated individuals [15]. Still, genomic copies/ml of the sample remain very high (~million genomic copies/ml), significantly above the LOD of this test, implicating that test will be useful in a post-vaccination era. We found some of the CLIA-defined LDT performance criteria difficult to apply. For example, CLIA would allow for up to 10% of base calls to be incorrect for accuracy determination, which, in the case of the ~30,000 bp SARS-CoV-2 genome, would mean ~3,000 inaccurate bases, which could lead to false variant detection. We accepted a minimum of 90% genome coverage and >200X median depth, but detection of the genomic variant was considered for final accuracy calculations. Because one erroneous SNP is unlikely to change genotyping conclusions in most instances, analysis was limited to overall variant detection using the default parameters for ease of implementation in the clinical laboratory. We also did not test the recommended 20 replicates to determine the LOD because this test is not intended to detect an analyte but the variations (genomic variant) in the analyte. Acceptable depth of coverage has been identified as 10X coverage of >90% of the genome. Low input of the RNA or lower reads will not meet these criteria. Therefore, false or undetermined variants are unlikely to be reported in low-input samples. Implementing a continuous performance measurement plan via an internal or external PT program is required to successfully integrate any test in the clinical laboratory (CAP Checklist 2021; <https://www.cap.org>). A set of reference SARS-CoV-2 variants is amenable to internal and external quality assurance testing. We assessed the entire workflow in preliminary internal PT by re-testing blind samples and inter-personnel reproducibility (details not shown). WGS is a dynamic technology evolving rapidly; therefore, our validated pipeline is unlikely to remain static. Re-validation provision is crucial for the seamless and timely implementation of changes to wet-lab reagents or the analysis pipeline. We have introduced a provision for reagent

verification at each lot change by re-testing samples in triplicate. Likewise, raw sequencing data will be re-analyzed with an updated analysis pipeline, and the accuracy of variant detection will be verified. DRAGEN COVID Lineage variant pipeline has been updated during the validation, and the variant identified by the two versions are in 100% concordance (data not shown). In general, WGS diagnosis reports are complex, and the format could pose challenges for the end-user. We adopted a simple format already used for SARS-CoV-2 diagnosis, and an updated report with the variant information will be issued if the reflex testing for variant detection is requested. This study possesses certain limitations. First, only a limited number of WGS-based assays were included in the validation study based on the limited application of this test for clinical decision-making. Second, we could not establish this test's clinical sensitivity and specificity because the clinical presentation of the patients infected with different variants is not distinct [30]. Moreover, the study could not acquire clinical samples of every lineage to demonstrate accuracy. However, CLIA and CAP regulations do not require validation of each mutation in the case of the mutation detection assay. For example, genomic mutation detection assays are commonly used in oncology. Likewise, we demonstrated the accuracy of the assay by testing three major variants. The vaccination status of the samples was not available to compare the application of the assay in the post-vaccination era. Although the vaccination status of the individual patient was not available, we tested >100 patient samples in the post-vaccination era to demonstrate that the test remains applicable to the vaccinated population. We could not try the test's clinical utility because that would require enrolling the patients infected with different variants and administering variant-specific treatment. Since their inception, most NGS-based testing has been limited to large medical centers, public health laboratories, or centralized genomics facilities with rather large infrastructures. The recent pandemic has accentuated the importance of de-centralized independent laboratories. For example, Advanta Genetics has served East Texas by testing > 500,000 SARS-CoV-2 samples. Thus, this validation can be used as guidelines for other small laboratories with NGS capacities if a need for SARS-CoV-2 variant detection arises. Although inevitable in the early stages, de-centralized NGS testing presents several challenges, such as high cost and turnaround time because of low volume testing. However, de-centralized and rapid testing for circulating SARS-CoV-2 variants may become crucial for clinical management and tracking the transmission at the local and regional levels. Sequencing cost in terms of dollars/gigabases has plummeted with high throughput instruments such as Illumina NovaSeq. However, such an instrument alone costs ~ a million dollars. It would require batching of ~30,000 samples to achieve the highest efficiency, which is not practical for independent laboratories, potentially leaving a gap in underserved communities.



This study introduces pre-pooling normalization to improve sequencing efficiency, which is crucial for smaller laboratories with low throughput sequencers. The emergence of more affordable sequencers such as Oxford Nanopore (Starting cost of \$10,000) enriches the opportunity for de-centralized genomic testing if a variant with distinct clinical needs emerges or for any future pandemic. We have demonstrated that NGS services, including clinical testing, could be delivered locally with well-defined quality metrics at an affordable cost. Global NGS data aggregators that emerged from this pandemic have been helpful for analysis support needed for resource-limited laboratories (<https://www.gisaid.org/collaborations/enabled-by-hcov-19-data-from-gisaid/>), but sequencing infrastructure remains centralized mainly [31]. The local-delivery model would also be more responsive to the target clients' needs and enhance the adoption of NGS across health care systems. We have demonstrated the application of this approach in the East Texas region and tracked the variant evolution throughout the pandemic. An alternate hybrid model has been proposed with complementary central and local services to balance the need for speed and investment [32]. The FDA Genome Tracker network for tracking foodborne pathogens and the Centers for Disease Control and Prevention (CDC) Advanced Molecular Detection (AMD) initiative for improving infectious disease surveillance are existing hybrid models in the United States [33, 34]. Notably, there are still significant challenges to implementing comprehensive WGS services locally [35, 36]. This study has established the performance specifications for NGS-based SARS-CoV-2 variant detection according to CAP and CLIA guidelines. We anticipate that the COVIDSeq LDT validation framework presented in this study, in synergy with increasingly accessible analysis support, will advance the localization of comprehensive NGS services in independent clinical laboratories. We have benchmarked quality assurance and quality control measures for implementing such testing and a simplified reporting format for end-users with limited NGS understanding. The study also affirmed the application of de-centralized NGS testing for clinical and public health applications with any resurgence of COVID-19 or the next infectious disease outbreak.

#### Author Contributions

RC: Conceived the study and assisted in the writing; VT: Performed the bioinformatics analysis and organized the data; SA: performed the experiments and organized the data; EB: performed the experiments; RS\*: Designed the study, supervised experiments and wrote the manuscript. RC, VT, and SA contributed equally to this work and shared the first authorship.

#### Ethical Statement

The study was exempted by IRB (Institutional Review Board) because only de-identified samples were used.

#### Data Availability

GISAID Identifier: EPI\_SET\_20220715vh doi: 10.55876/gis8.220715vh

All genome sequences and associated metadata in this dataset are published in GISAID's EpiCoV database. To view the contributors of each sequence with details such as accession number, Virus name, Collection date, Originating Lab and Submitting Lab, and the list of Authors, visit [10.55876/gis8.220715vh](https://gisaid.org/EPI_SET_20220715vh)

#### Data Snapshot

EPI\_SET\_20220715vh is composed of 72 individual genome sequences. The collection dates range from 2020-08-01 to 2022-09-25; Data were collected in 1 country and territory; All sequences in this dataset are compared to hCoV-19/Wuhan/WIV04/2019 (WIV04), the official reference sequence employed by GISAID (EPI\_ISL\_402124). Learn more at <https://gisaid.org/WIV04>.

#### References

1. Chilamakuri R, Agarwal S. COVID-19: Characteristics and Therapeutics. *Cells* 10 (2021): 206.
2. Lekana-Douki SE, N'dilimabaka N, Levasseur A, et al. Screening and Whole Genome Sequencing of SARS-CoV-2 Circulating During the First Three Waves of the COVID-19 Pandemic in Libreville and the Haut-Ogooué Province in Gabon. *Front Med (Lausanne)* 9 (2022): 877391.
3. John G, Sahajpal NS, Mondal AK, et al. Next-Generation Sequencing (NGS) in COVID-19: A Tool for SARS-CoV-2 Diagnosis, Monitoring New Strains and Phylodynamic Modeling in Molecular Epidemiology. *Curr Issues Mol Biol* 43 (2021): 845-867.
4. Galloway SE, Paul P, MacCannell DR, et al. Emergence of SARS-CoV-2 B.1.1.7 Lineage - United States, December 29, 2020-January 12, 2021. *MMWR Morb Mortal Wkly Rep* 70 (2021): 95-99.
5. Cascella M, Rajnik M, Aleem A, et al. Features, Evaluation, and Treatment of Coronavirus (COVID-19). 2022 Jun 30. In: StatPearls [Internet]. Treasure Island (FL): StatPearls Publishing (2022).
6. Ibn-Mohammed T, Mustapha KB, Godsell J, et al. A critical analysis of the impacts of COVID-19 on the global economy and ecosystems and opportunities for circular economy strategies. *Resour Conserv Recycl* 164 (2021): 105169.
7. Banoun H. Evolution of SARS-CoV-2: Review of Mutations, Role of the Host Immune System. *Nephron* 145 (2021): 392-403.
8. Aleem A, Akbar Samad AB, Slenker AK. Emerging

- Variants of SARS-CoV-2 and Novel Therapeutics Against Coronavirus (COVID-19). 2022 May 12. In: StatPearls [Internet]. Treasure Island (FL): StatPearls Publishing (2022).
9. Iketani S, Liu L, Guo Y, et al. Antibody evasion properties of SARS-CoV-2 Omicron sublineages. *Nature* 604 (2022): 553-556.
  10. Jo DH, Minn D, Lim J, et al. Rapidly Declining SARS-CoV-2 Antibody Titers within 4 Months after BNT162b2 Vaccination. *Vaccines (Basel)* 9 (2021): 1145.
  11. Levin EG, Lustig Y, Cohen C, et al. Waning Immune Humoral Response to BNT162b2 Covid-19 Vaccine over 6 Months. *N Engl J Med* 385 (2021): e84.
  12. Chen RE, Zhang X, Case JB, et al. Resistance of SARS-CoV-2 variants to neutralization by monoclonal and serum-derived polyclonal antibodies. *Nat Med* 27 (2021): 717-726.
  13. Shen X, Tang H, Pajon R, et al. Neutralization of SARS-CoV-2 Variants B.1.429 and B.1.351. *N Engl J Med* 384 (2021): 2352-2354.
  14. Carpenter RE, Tamrakar V, Chahar H, et al. Confirming Multiplex RT-qPCR Use in COVID-19 with Next-Generation Sequencing: Strategies for Epidemiological Advantage. *Glob Health Epidemiol Genom* (2022): 2270965.
  15. Jennings LJ, Arcila ME, Corless C, et al. Guidelines for Validation of Next-Generation Sequencing-Based Oncology Panels: A Joint Consensus Recommendation of the Association for Molecular Pathology and College of American Pathologists. *J Mol Diagn* 19 (2017): 341-365.
  16. Puhach O, Adea K, Hulo N, et al. Infectious viral load in unvaccinated and vaccinated individuals infected with ancestral, Delta or Omicron SARS-CoV-2. *Nat Med* 28 (2022): 1491-1500.
  17. Khare S, Gurry C, Freitas L, et al. GISAID's Role in Pandemic Response. *China CDC Wkly* 3 (2021): 1049-1051.
  18. Charre C, Ginevra C, Sabatier M, et al. Evaluation of NGS-based approaches for SARS-CoV-2 whole genome characterisation. *Virus Evol* 6 (2020): veaa075.
  19. Welch NL, Zhu M, Hua C, et al. Multiplexed CRISPR-based microfluidic platform for clinical testing of respiratory viruses and identification of SARS-CoV-2 variants. *Nat Med* 28 (2022): 1083-1094.
  20. Pillay S, Giandhari J, Tegally H, et al. Whole Genome Sequencing of SARS-CoV-2: Adapting Illumina Protocols for Quick and Accurate Outbreak Investigation during a Pandemic. *Genes (Basel)* 11 (2020): 949.
  21. Bhoyar RC, Jain A, Sehgal P, et al. High throughput detection and genetic epidemiology of SARS-CoV-2 using COVIDSeq next-generation sequencing. *PLoS One* 16 (2021): e0247115.
  22. Huai Luo C, Paul Morris C, Sachithanandham J, et al. Infection With the Severe Acute Respiratory Syndrome Coronavirus 2 (SARS-CoV-2) Delta Variant Is Associated With Higher Recovery of Infectious Virus Compared to the Alpha Variant in Both Unvaccinated and Vaccinated Individuals. *Clin Infect Dis* 75 (2022): e715-e725.
  23. Cele S, Jackson L, Khoury DS, et al. Omicron extensively but incompletely escapes Pfizer BNT162b2 neutralization. *Nature* 602 (2022): 654-656.
  24. Gargis AS, Kalman L, Berry MW, et al. Assuring the quality of next-generation sequencing in clinical laboratory practice. *Nat Biotechnol* 30 (2012): 1033-1036.
  25. Kalman LV, Datta V, Williams M, et al. Development and Characterization of Reference Materials for Genetic Testing: Focus on Public Partnerships. *Ann Lab Med* 36 (2016): 513-520.
  26. Ou J, Lan W, Wu X, et al. Tracking SARS-CoV-2 Omicron diverse spike gene mutations identifies multiple inter-variant recombination events. *Signal Transduct Target Ther* 7 (2022): 138.
  27. Maxmen A. One million coronavirus sequences: popular genome site hits mega milestone. *Nature* 593 (2021): 21.
  28. Lindsey RL, Pouseele H, Chen JC, et al. Implementation of Whole Genome Sequencing (WGS) for Identification and Characterization of Shiga Toxin-Producing *Escherichia coli* (STEC) in the United States. *Front Microbiol* 7 (2016): 766.
  29. Walker TM, Kohl TA, Omar SV, et al. Whole-genome sequencing for prediction of *Mycobacterium tuberculosis* drug susceptibility and resistance: a retrospective cohort study. *Lancet Infect Dis* 15 (2015): 1193-1202.
  30. Hoang VT, Colson P, Levasseur A, et al. Clinical outcomes in patients infected with different SARS-CoV-2 variants at one hospital during three phases of the COVID-19 epidemic in Marseille, France. *Infect Genet Evol* 95 (2021): 105092.
  31. Udugama B, Kadhireshan P, Kozlowski HN, et al. Diagnosing COVID-19: The Disease and Tools for Detection. *ACS Nano* 14 (2020): 3822-3835.
  32. Arnold C. Considerations in centralizing whole genome sequencing for microbiology in a public health setting. *Expert Rev Mol Diagn* 16 (2016): 619-621.

33. Auffray C, Caulfield T, Griffin JL, et al. From genomic medicine to precision medicine: highlights of 2015. *Genome Med* 8 (2016): 12.
34. Allard MW, Strain E, Melka D, et al. Practical Value of Food Pathogen Traceability through Building a Whole-Genome Sequencing Network and Database. *J Clin Microbiol* 54 (2016): 1975-1983.
35. Lesho E, Clifford R, Onmus-Leone F, et al. The Challenges of Implementing Next Generation Sequencing Across a Large Healthcare System, and the Molecular Epidemiology and Antibiotic Susceptibilities of Carbapenemase-Producing Bacteria in the Healthcare System of the U.S. Department of Defense. *PLoS One* 11 (2016): e 0155770.
36. Robilotti E, Kamboj M. Integration of whole-genome sequencing into infection control practices: the potential and the hurdles. *J Clin Microbiol* 53 (2015): 1054-1055.

**Supplementary Table 1:** Demographic information of all 162 samples of SARS-COV-2, East Texas Region.

Accession ID	GISAID Accession ID	Date of Collection	County	State	Run Date	PANGO Lineage
ACSQ1-12	EPI_ISL_8719355	08/01/2020	Smith	TX	9/15/21	B.1
ACSQ1-8	EPI_ISL_8729880	08/01/2020	Smith	TX	9/15/21	B.1.243
ACSQ1-13	EPI_ISL_8729884	08/01/2020	Smith	TX	9/15/21	B.1.234
ACSQ1-14	EPI_ISL_8729885	08/01/2020	Smith	TX	9/15/21	B.1.126
ACSQ1-15	EPI_ISL_8729886	08/01/2020	Smith	TX	9/15/21	B.1.602
ACSQ1-16	EPI_ISL_8729887	08/01/2020	Smith	TX	9/15/21	B.1
ACSQ1-17	EPI_ISL_8729888	08/01/2020	Smith	TX	9/15/21	B.1
ACSQ1-18	EPI_ISL_8729889	08/01/2020	Smith	TX	9/15/21	B.1.564
ACSQ1-19	EPI_ISL_8729890	08/01/2020	Smith	TX	9/15/21	B.1.617.2
ACSQ1-20	EPI_ISL_8729891	08/01/2020	Smith	TX	9/15/21	AY.3
ACSQ1-9	EPI_ISL_8729881	08/01/2020	Smith	TX	9/15/21	B.1.574
ACSQ1-21	EPI_ISL_8729892	08/01/2020	Smith	TX	9/15/21	B.1.617.2
ACSQ1-10	EPI_ISL_8729882	08/01/2020	Smith	TX	9/15/21	B.1.574
ACSQ1-11	EPI_ISL_8729883	08/01/2020	Smith	TX	9/15/21	B.1.2
ACSQ1-3	EPI_ISL_8729875	7/27/21	Smith	TX	9/15/21	AY.3
FT-SA99822	EPI_ISL_3523224	7/27/21	Smith	TX	08/11/2021	AY.3
FT-SA99824	EPI_ISL_3523226	7/27/21	Smith	TX	08/11/2021	AY.3
FT-SA99827	EPI_ISL_3523229	7/27/21	Smith	TX	08/11/2021	AY.3
FT-SA99828	EPI_ISL_3523230	7/27/21	Brazos	TX	08/11/2021	B.1.617.2
FT-SA99833	EPI_ISL_3523235	7/27/21	Smith	TX	08/11/2021	AY.3
ACSQ1-5	EPI_ISL_8729877	7/28/21	Smith	TX	9/15/21	AY.3
FT-SA99772	EPI_ISL_3523183	7/28/21	Tarrant	TX	08/11/2021	B.1.617.2
FT-SA99802	EPI_ISL_3523208	7/28/21	Van zandt	TX	08/11/2021	AY.3
FT-SA99820	EPI_ISL_3523222	7/28/21	Smith	TX	08/11/2021	AY.3
FT-SA99821	EPI_ISL_3523223	7/28/21	Smith	TX	08/11/2021	AY.3
FT-SA99823	EPI_ISL_3523225	7/28/21	Smith	TX	08/11/2021	AY.3
FT-SA99825	EPI_ISL_3523227	7/28/21	Smith	TX	08/11/2021	AY.3
FT-SA99826	EPI_ISL_3523228	7/28/21	Smith	TX	08/11/2021	AY.3
FT-SA99829	EPI_ISL_3523231	7/28/21	Smith	TX	08/11/2021	AY.3
FT-SA99830	EPI_ISL_3523232	7/28/21	Smith	TX	08/11/2021	B.1.617.2
FT-SA99832	EPI_ISL_3523234	7/28/21	Johnson	TX	08/11/2021	AY.3
FT-SA99852	EPI_ISL_3523252	7/28/21	Smith	TX	08/11/2021	AY.3
FT-SA99856	EPI_ISL_3523256	7/28/21	Smith	TX	08/11/2021	B.1.617.2
FT-SA99861	EPI_ISL_3523261	7/28/21	Smith	TX	08/11/2021	B.1.617.2
FT-SA99863	EPI_ISL_3523262	7/28/21	Gregg	TX	08/11/2021	B.1.617.2
ACSQ1-7	EPI_ISL_8729879	7/29/21	Smith	TX	9/15/21	AY.3
ACSQ1-2	EPI_ISL_8729874	7/29/21	Smith	TX	9/15/21	AY.25
FT-SA99774	EPI_ISL_3523185	7/29/21	Hidalgo	TX	08/11/2021	B.1.617.2
FT-SA99785	EPI_ISL_3523194	7/29/21	Smith	TX	08/11/2021	AY.3
FT-SA99790	EPI_ISL_3523198	7/29/21	Smith	TX	08/11/2021	AY.3
FT-SA99791	EPI_ISL_3523199	7/29/21	Smith	TX	08/11/2021	AY.3
FT-SA99792	EPI_ISL_3523200	7/29/21	Harris	TX	08/11/2021	B.1.617.2
FT-SA99793	EPI_ISL_3523201	7/29/21	Smith	TX	08/11/2021	AY.3
FT-SA99795	EPI_ISL_3523202	7/29/21	Smith	TX	08/11/2021	AY.3



FT-SA99796	EPI_ISL_3523203	7/29/21	Johnson	TX	08/11/2021	B.1.617.2
FT-SA99803	EPI_ISL_3523209	7/29/21	Smith	TX	08/11/2021	AY.3
FT-SA99804	EPI_ISL_3523210	7/29/21	Hood	TX	08/11/2021	B.1.617.2
FT-SA99831	EPI_ISL_3523233	7/29/21	Tarrant	TX	08/11/2021	B.1.617.2
FT-SA99835	EPI_ISL_3523237	7/29/21	Smith	TX	08/11/2021	AY.3
FT-SA99840	EPI_ISL_3523241	7/29/21	Smith	TX	08/11/2021	AY.3
FT-SA99841	EPI_ISL_3523242	7/29/21	Smith	TX	08/11/2021	AY.3
FT-SA99850	EPI_ISL_3523250	7/29/21	Smith	TX	08/11/2021	AY.3
FT-SA99858	EPI_ISL_3523258	7/29/21	Smith	TX	08/11/2021	AY.3
FT-SA99859	EPI_ISL_3523259	7/29/21	Smith	TX	08/11/2021	AY.3
FT-SA99860	EPI_ISL_3523260	7/29/21	Smith	TX	08/11/2021	B.1.617.2
FT-SA99865	EPI_ISL_3523264	7/29/21	Smith	TX	08/11/2021	AY.3
ACSQ1-4	EPI_ISL_8729876	7/30/21	Smith	TX	9/15/21	AY.3
FT-SA99770	EPI_ISL_3523181	7/30/21	Smith	TX	08/11/2021	B.1.617.2
FT-SA99771	EPI_ISL_3523182	7/30/21	Smith	TX	08/11/2021	AY.3
FT-SA99773	EPI_ISL_3523184	7/30/21	Brazos	TX	08/11/2021	B.1.617.2
FT-SA99775	EPI_ISL_3523186	7/30/21	Smith	TX	08/11/2021	AY.3
FT-SA99776	EPI_ISL_3523187	7/30/21	Harrison	TX	08/11/2021	AY.3
FT-SA99777	EPI_ISL_3523188	7/30/21	Smith	TX	08/11/2021	AY.3
FT-SA99780	EPI_ISL_3523189	7/30/21	Smith	TX	08/11/2021	B.1.617.2
FT-SA99781	EPI_ISL_3523190	7/30/21	Smith	TX	08/11/2021	B.1.617.2
FT-SA99782	EPI_ISL_3523191	7/30/21	Smith	TX	08/11/2021	AY.3
FT-SA99783	EPI_ISL_3523192	7/30/21	Smith	TX	08/11/2021	AY.3
FT-SA99784	EPI_ISL_3523193	7/30/21	Smith	TX	08/11/2021	AY.3
FT-SA99786	EPI_ISL_3523195	7/30/21	Smith	TX	08/11/2021	AY.3
FT-SA99788	EPI_ISL_3523196	7/30/21	Smith	TX	08/11/2021	AY.3
FT-SA99789	EPI_ISL_3523197	7/30/21	Smith	TX	08/11/2021	AY.3
FT-SA99797	EPI_ISL_3523204	7/30/21	Smith	TX	08/11/2021	AY.3
FT-SA99799	EPI_ISL_3523205	7/30/21	Smith	TX	08/11/2021	AY.3
FT-SA99800	EPI_ISL_3523206	7/30/21	Smith	TX	08/11/2021	B.1.617.2
FT-SA99801	EPI_ISL_3523207	7/30/21	Henderson	TX	08/11/2021	B.1.617.2
FT-SA99805	EPI_ISL_3523211	7/30/21	Smith	TX	08/11/2021	B.1.617.2
FT-SA99809	EPI_ISL_3523214	7/30/21	Smith	TX	08/11/2021	B.1.617.2
ACSQ1-6	EPI_ISL_8729878	7/31/21	Smith	TX	9/15/21	B.1.617.2
FT-SA99807	EPI_ISL_3523212	7/31/21	Smith	TX	08/11/2021	B.1.617.2
FT-SA99808	EPI_ISL_3523213	7/31/21	Smith	TX	08/11/2021	B.1.617.2
FT-SA99813	EPI_ISL_3523215	7/31/21	Johnson	TX	08/11/2021	B.1.617.2
FT-SA99815	EPI_ISL_3523217	7/31/21	Smith	TX	08/11/2021	AY.3
FT-SA99816	EPI_ISL_3523218	7/31/21	Smith	TX	08/11/2021	AY.3
FT-SA99817	EPI_ISL_3523219	7/31/21	Smith	TX	08/11/2021	AY.3
FT-SA99818	EPI_ISL_3523220	7/31/21	Smith	TX	08/11/2021	AY.3
FT-SA99819	EPI_ISL_3523221	7/31/21	Henderson	TX	08/11/2021	B.1.617.2
FT-SA99814	EPI_ISL_3523216	08/11/2021	Smith	TX	08/11/2021	AY.3
ACSQ4-23	EPI_ISL_8629675	12/16/21	Passaic	NJ	12/26/21	B.1.1.529
ACSQ4-22	EPI_ISL_8629674	12/17/21	Smith	TX	12/26/21	B.1.1.529
ACSQ4-21	EPI_ISL_8629673	12/17/21	Tarrant	TX	12/26/21	AY.111

ACSQ4-19	EPI_ISL_8629671	12/18/21	Smith	TX	12/26/21	BA.1
ACSQ4-1	EPI_ISL_8629666	12/20/21	Tarrant	TX	12/26/21	B.1.1.529
ACSQ4-58	EPI_ISL_8428058	12/20/21	Smith	TX	12/26/21	BA.1
ACSQ4-4	EPI_ISL_8629668	12/20/21	Smith	TX	12/26/21	AY.103
ACSQ4-18	EPI_ISL_8629670	12/21/21	Navarro	TX	12/26/21	B.1.1.529
ACSQ4-2	EPI_ISL_8629667	12/21/21	Navarro	TX	12/26/21	B.1.1.529
ACSQ4-20	EPI_ISL_8629672	12/22/21	Navarro	TX	12/26/21	AY.3
ACSQ4-24	EPI_ISL_8629676	12/22/21	Passaic	NJ	12/26/21	B.1.1.529
ACSQ4-17	EPI_ISL_8629669	12/22/21	Smith	TX	12/26/21	B.1.1.529
ACSQ7-1	EPI_ISL_13655593	4/18/22	Van zandt	TX	6/24/22	BA.2
ACSQ7-2	EPI_ISL_13655594	4/19/22	Smith	TX	6/24/22	BA.3
ACSQ7-11	EPI_ISL_13655603	4/28/22	Smith	TX	6/24/22	BA.2.3
ACSQ7-17	EPI_ISL_13655607	05/01/2022	Brazos	TX	6/24/22	BA.5
ACSQ7-15	EPI_ISL_13671513	05/09/2022	Brazos	TX	6/24/22	BA.2.9
ACSQ7-3	EPI_ISL_13655595	05/10/2022	Smith	TX	6/24/22	BA.2.12.1
ACSQ7-12	EPI_ISL_13655604	05/10/2022	Brazos	TX	6/24/22	BA.2.3
ACSQ7-16	EPI_ISL_13655606	05/11/2022	Van zandt	TX	6/24/22	BA.2.9
ACSQ7-5	EPI_ISL_13655597	5/16/22	Smith	TX	6/24/22	BA.2.12.1
ACSQ7-4	EPI_ISL_13655596	5/16/22	Smith	TX	6/24/22	BA.2.12.1
ACSQ7-7	EPI_ISL_13655599	5/17/22	Tarrant	TX	6/24/22	BA.2.12.1
ACSQ7-6	EPI_ISL_13655598	5/18/22	Smith	TX	6/24/22	BA.2.12.1
ACSQ8-1	EPI_ISL_13762660	5/18/22	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ7-13	EPI_ISL_13671512	5/20/22	Smith	TX	6/24/22	BA.2.37
ACSQ7-8	EPI_ISL_13655600	5/21/22	Henderson	TX	6/24/22	BA.2.12.1
ACSQ7-14	EPI_ISL_13655605	5/23/22	Johnson	TX	6/24/22	BA.2.37
ACSQ8-2	EPI_ISL_13819250	06/02/2022	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ7-10	EPI_ISL_13655602	06/02/2022	Smith	TX	6/24/22	BA.2.12.1
ACSQ7-9	EPI_ISL_13655601	06/02/2022	Smith	TX	6/24/22	BA.2.12.1
ACSQ8-3	EPI_ISL_13819251	6/21/22	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-4	EPI_ISL_13762661	6/21/22	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-5	EPI_ISL_13762662	6/24/22	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-7	EPI_ISL_13762664	6/24/22	Gregg	TX	07/09/2022	BA.5.2.1
ACSQ8-6	EPI_ISL_13762663	6/25/22	Gregg	TX	07/09/2022	BA.5.5
ACSQ8-8	EPI_ISL_13762665	6/25/22	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-9	EPI_ISL_13762666	6/26/22	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-10	EPI_ISL_13762667	6/26/22	Gregg	TX	07/09/2022	BA.4.1
ACSQ8-11	EPI_ISL_13762668	6/28/22	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-12	EPI_ISL_13819252	6/30/22	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-14	EPI_ISL_13819253	6/30/22	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-13	EPI_ISL_13762669	07/01/2022	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-15	EPI_ISL_13762670	07/01/2022	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-16	EPI_ISL_13762671	07/01/2022	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-18	EPI_ISL_13762673	07/03/2022	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-17	EPI_ISL_13762672	07/04/2022	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-19	EPI_ISL_13762674	07/05/2022	Gregg	TX	07/09/2022	BA.5.2.1
ACSQ8-22	EPI_ISL_13762677	07/05/2022	Gregg	TX	07/09/2022	BA.2.12.1

ACSQ8-23	EPI_ISL_13762678	07/05/2022	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ8-20	EPI_ISL_13762675	07/06/2022	Gregg	TX	07/09/2022	BA.5.2.1
ACSQ8-21	EPI_ISL_13762676	07/06/2022	Gregg	TX	07/09/2022	BA.2.12.1
ACSQ9-15	EPI_ISL_15169383	09/10/2022	Smith	TX	9/26/22	BA.5.2
ACSQ9-24	EPI_ISL_15169391	09/10/2022	Henderson	TX	9/26/22	BA.5.2.1
ACSQ9-23	EPI_ISL_15169390	09/11/2022	Smith	TX	9/26/22	BA.2.3
ACSQ9-22	EPI_ISL_15169389	09/11/2022	Anderson	TX	9/26/22	BA.5.5
ACSQ9-17	EPI_ISL_15169385	09/12/2022	Tarrant	TX	9/26/22	BA.5.6
ACSQ9-18	EPI_ISL_15169386	9/13/22	Henderson	TX	9/26/22	BA.5.1.6
ACSQ9-21	EPI_ISL_15169388	9/13/22	Shelby	TX	9/26/22	BA.5.2.1
ACSQ9-14	EPI_ISL_15169382	9/13/22	Tarrant	TX	9/26/22	BA.5.2
ACSQ9-20	EPI_ISL_15169387	9/14/22	Nacogdoches	TX	9/26/22	BA.5.2
ACSQ9-13	EPI_ISL_15169381	9/14/22	Smith	TX	9/26/22	BA.4.6
ACSQ9-12	EPI_ISL_15169380	9/14/22	Smith	TX	9/26/22	BA.4.6
ACSQ9-7	EPI_ISL_15169375	9/16/22	Smith	TX	9/26/22	BE.1.1
ACSQ9-16	EPI_ISL_15169384	9/16/22	Anderson	TX	9/26/22	BA.5.5
ACSQ9-8A	EPI_ISL_15169376	9/16/22	Denton	TX	9/26/22	BA.5.2
ACSQ9-9	EPI_ISL_15169377	9/17/22	Smith	TX	9/26/22	BE.1.1
ACSQ9-11	EPI_ISL_15169379	9/18/22	Brazos	TX	9/26/22	BA.4.6
ACSQ9-10	EPI_ISL_15169378	9/18/22	Gregg	TX	9/26/22	BA.5.2
ACSQ9-3	EPI_ISL_15169371	9/19/22	Johnson	TX	9/26/22	BA.5.1.1
ACSQ9-5	EPI_ISL_15169373	9/19/22	Tarrant	TX	9/26/22	BA.5.1
ACSQ9-4A	EPI_ISL_15169372	9/19/22	Tarrant	TX	9/26/22	BA.5.2.1
ACSQ9-2	EPI_ISL_15169370	9/20/22	Smith	TX	9/26/22	BA.5.2
ACSQ9-6	EPI_ISL_15169374	9/21/22	Henderson	TX	9/26/22	BA.5.1
ACSQ9-1A	EPI_ISL_15169369	9/22/22	Smith	TX	9/26/22	BA.4.6

**Supplementary Table 2:** Intra-Sequencing run precision: Variant identification concordance among the triplicate testing of the same samples.

\*Variant detection was inconclusive because of the low coverage.

Accession	(Rep-1)				(Rep-2)				(Rep-3)				Concordance
	Median Coverage	Coverage >= 30x	Pango Lineage	WHO label	Median Coverage	Coverage >= 30x	Pango Lineage	WHO label	Median Coverage	Coverage >= 30x	Pango Lineage	WHO label	
AACSQ1-3	345	86.06%	AY.3	Delta (B.1.617.2-like)	302	86.55%	AY.3	Delta (B.1.617.2-like)	309	85.61%	AY.3	Delta (B.1.617.2-like)	YES
AACSQ1-3*	120	80.09%	AY.3	Delta (B.1.617.2-like)	108	81.39%	None		107	79.47%	AY.3	Delta (B.1.617.2-like)	Inconclusive
AACSQ1-3*	149	84.77%	AY.82	Delta (B.1.617.2-like)	132	84.38%	None		85	74.01%	AY.3	Delta (B.1.617.2-like)	Inconclusive
ACSQ1-7	2289	99.78%	AY.25	Delta (B.1.617.2-like)	684	99.19%	AY.25	Delta (B.1.617.2-like)	773	99.39%	AY.25	Delta (B.1.617.2-like)	YES
ACSQ1-7	364	99.33%	AY.25	Delta (B.1.617.2-like)	289	99.14%	AY.25	Delta (B.1.617.2-like)	371	98.61%	AY.25	Delta (B.1.617.2-like)	YES
ACSQ1-7	444	98.70%	AY.25	Delta (B.1.617.2-like)	303	98.43%	AY.25	Delta (B.1.617.2-like)	435	98.64%	AY.25	Delta (B.1.617.2-like)	YES
ACSQ1-2	338	98.92%	AY.3	Delta (B.1.617.2-like)	349	98.97%	AY.3	Delta (B.1.617.2-like)	316	98.92%	AY.3	Delta (B.1.617.2-like)	YES
ACSQ1-2	418	98.34%	AY.3	Delta (B.1.617.2-like)	427	98.45%	AY.3	Delta (B.1.617.2-like)	336	98.30%	AY.3	Delta (B.1.617.2-like)	YES
ACSQ2-14	547	92.60%	AY.3	Delta (B.1.617.2-like)	446	92.25%	AY.3	Delta (B.1.617.2-like)	560	93.01%	AY.3	Delta (B.1.617.2-like)	YES
ACSQ2-14	226	92.85%	AY.3	Delta (B.1.617.2-like)	165	91.83%	AY.3	Delta (B.1.617.2-like)	224	93.91%	AY.3	Delta (B.1.617.2-like)	YES
ACSQ2-14	241	94.21%	AY.3	Delta (B.1.617.2-like)	227	94.18%	AY.3	Delta (B.1.617.2-like)	242	94.38%	AY.3	Delta (B.1.617.2-like)	YES
ACSQ2-20	757	99.59%	AY.39.1	Delta (B.1.617.2-like)	706	99.59%	AY.39.1	Delta (B.1.617.2-like)	786	99.59%	AY.39.1	Delta (B.1.617.2-like)	YES
ACSQ2-20	310	99.39%	AY.39.1	Delta (B.1.617.2-like)	295	99.37%	AY.39.1	Delta (B.1.617.2-like)	330	99.70%	AY.39.1	Delta (B.1.617.2-like)	YES
ACSQ2-20	383	99.01%	AY.39.1	Delta (B.1.617.2-like)	371	98.89%	AY.39.1	Delta (B.1.617.2-like)	410	98.92%	AY.39.1	Delta (B.1.617.2-like)	YES
ACSQ2-17	633	95.39%	AY.103	Delta (B.1.617.2-like)	583	95.43%	AY.103	Delta (B.1.617.2-like)	595	95.20%	AY.103	Delta (B.1.617.2-like)	YES
ACSQ2-17	265	95.65%	AY.103	Delta (B.1.617.2-like)	237	94.94%	AY.103	Delta (B.1.617.2-like)	263	95.17%	AY.103	Delta (B.1.617.2-like)	YES
ACSQ2-17	298	97.21%	AY.103	Delta (B.1.617.2-like)	275	96.96%	AY.103	Delta (B.1.617.2-like)	247	96.38%	AY.103	Delta (B.1.617.2-like)	YES
ACSQ2-AB	1781	98.24%	BA.5.2	Omicron (BA.5-like)	1249	96.34%	BA.5.2	Omicron (BA.5-like)	2050	97.95%	BA.5.2	Omicron (BA.5-like)	YES



**Supplementary Table 3:** Inter-Sequencing run precision: Variant identification concordance across three different library preparation and sequencing runs.

Accession	ACSQ-2 (Rep-1)				ACSQ-3 (Rep-2)				ACSQ-4 (Rep-3)				Concordance
	Median Coverage-R1	Coverage >= 30x-R1	Pango Lineage	WHO label	Median Coverage-R2	Coverage >= 30xR-2	Pango Lineage	WHO label	Median Coverage-3	Coverage >= 30x-R3	Pango Lineage	WHO label	
ACSQ1-3	345	86.06%	AY.3	Delta (B.1.617.2-like)	120	80.09%	AY.3	Delta (B.1.617.2-like)	149	84.77%	AY.82	Delta (B.1.617.2-like)	YES
ACSQ1-7	2289	99.78%	AY.25	Delta (B.1.617.2-like)	364	99.33%	AY.25	Delta (B.1.617.2-like)	444	98.70%	AY.25	Delta (B.1.617.2-like)	YES
ACSQ1-2	797	99.01%	AY.3	Delta (B.1.617.2-like)	338	98.92%	AY.3	Delta (B.1.617.2-like)	418	98.34%	AY.3	Delta (B.1.617.2-like)	YES
ACSQ2-14	547	92.60%	AY.3	Delta (B.1.617.2-like)	226	92.85%	AY.3	Delta (B.1.617.2-like)	241	94.21%	AY.3	Delta (B.1.617.2-like)	YES
ACSQ2-20	757	99.59%	AY.39.1	Delta (B.1.617.2-like)	310	99.39%	AY.39.1	Delta (B.1.617.2-like)	383	99.01%	AY.39.1	Delta (B.1.617.2-like)	YES
ACSQ2-17	633	95.39%	AY.103	Delta (B.1.617.2-like)	265	95.65%	AY.103	Delta (B.1.617.2-like)	298	97.21%	AY.103	Delta (B.1.617.2-like)	YES
ACSQ4-1	1286	98.32%	B.1.1.529	Omicron (B.1.1.529-like)	1363	99.15%	B.1.1.529	Omicron (B.1.1.529-like)	1209	97.54%	B.1.1.529	Omicron (B.1.1.529-like)	YES
ACSQ4-85	1131	96.64%	B.1.1.529	Omicron (B.1.1.529-like)	1201	97.10%	B.1.1.529	Omicron (B.1.1.529-like)	1308	96.61%	B.1.1.529	Omicron (B.1.1.529-like)	YES
ACSQ4-2	867	92.42%	B.1.1.529	Omicron (B.1.1.529-like)	1130	94.96%	B.1.1.529	Omicron (B.1.1.529-like)	900	95.46%	B.1.1.529	Omicron (B.1.1.529-like)	YES